


Ensemble Keyboard Arm

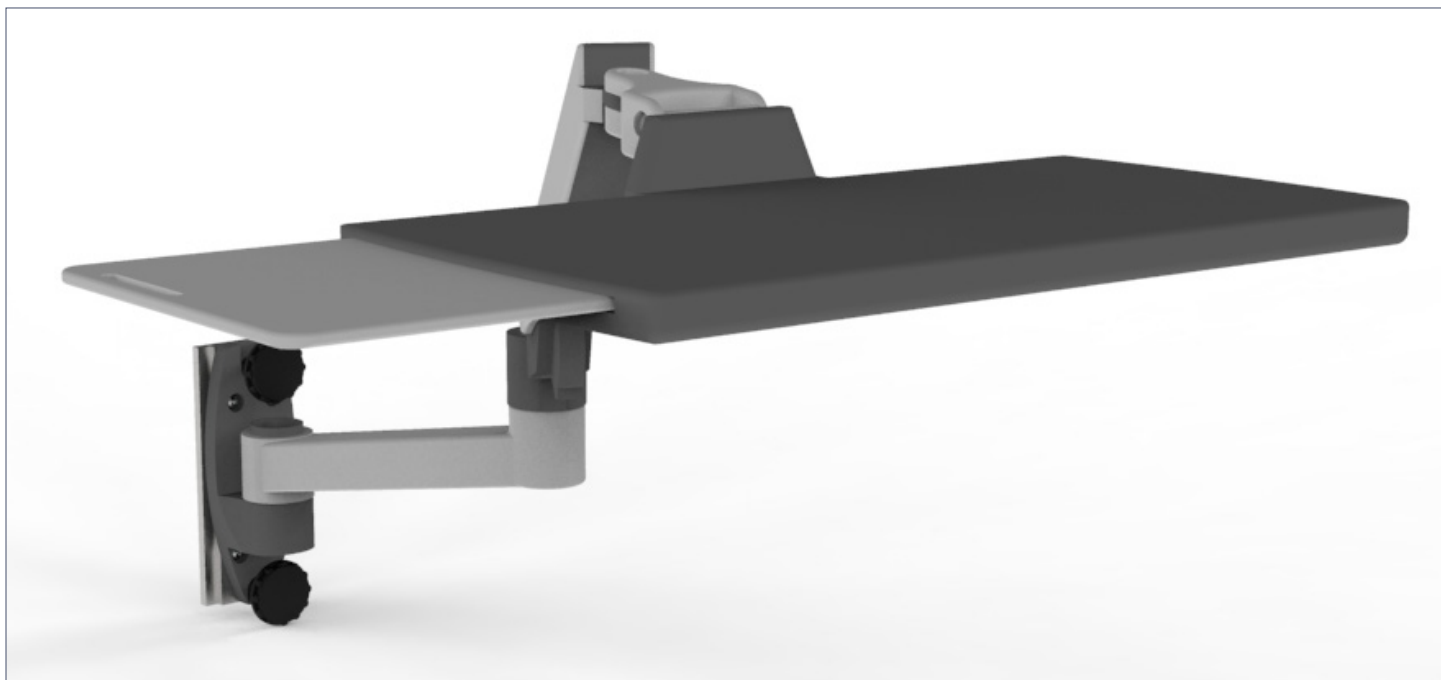

SPECIFICATIONS

Construction

- Designed to mount into the existing track on the column opposite the sink
- Once installed, track allows arm to be moved up and down the track by loosening two large, knurled knobs, without needing any tools

Capacity

- 5 lbs

Materials

- Keyboard arm components are powder-coated aluminum and ABS plastic for the keyboard tray

Dimensions

- 22.3" L x 33.5" D x 16.1" H
- Keyboard extended: 22.3" L (shown in photo)
- Keyboard not extended: 17.9" L

Performance

- Has several pivot points to articulate up and down, and side to side.
- Can be pulled out when in use, then pushed back when not in use
- Arm has 33" and 180° of sweep, and 13" of lift

Compatibility

- Fits all Mopec Maestro and Maestro Encore grossing stations

NOTE:

- If ordering with AX034 Monitor Arm, please order AX070 bundle part number.

MOPEC
GROUP

mopec
ELEVATING PATHOLOGY

MORTECH
MANUFACTURING

www.mopecgroup.com

800.362.8491

800 Tech Row, Madison Heights, MI